

Oxide Line - Valve Controller

Software Functional Requirements

By **Dr Gary Twinn**

Contents

Oxide Line - Valve Controller	1
Contents	2
List of Figures	3
List of Tables	4
Overview	5
Valve Controller Description	7
Input Messages	8
Output Messages	9
Valve to GPIO assignment	11
Driver circuits and power supply	12
Web interface	14
Interfaces	15
9 Pin Valve Cable Specification	15
9 Pin RS232 Cable specification	16
4 Pin RS485 Cable Specification	17
USB Connections.....	18
Gauge Data Request Commands.....	18
Raspberry Pi Operating System Installation.....	19
Operating system.....	19
Raspberry Pi Software installation	19
Run Operating System Updates.....	19
Install git and gh.....	19
Configure git and download the valve controller code	20
Create a python venv	20
Change to the venv.....	20
Install the Required Python Libraries.....	20
Gunicorn for Python service	20
Nginx installation.....	21

List of Figures

Figure 1: The UCL Oxide Nobel Gass Extraction Line.....	5
Figure 2: Schematic for a valve driver board	13
Figure 3: Web status page	14
Figure 4: Valve controller 9 pin D plug assignments.....	15
Figure 5: RS232 D Plug Assignment and null modem cable	16
Figure 6: RS485 M12 connector layout and intercept cable	17
Figure 7: USB connections to the Valve Controller Raspberry Pi	18

List of Tables

Table 1: Manual Valves.....	5
Table 2: Automated valves.....	6
Table 3: GPIO to valve assignments	11
Table 4: External Equipment Connections	18
Table 5: Commands to Vacuum Gauges.....	18

Overview

The new oxide extraction line is a configurable Nobel gas extraction line used as a test environment (Figure 1).

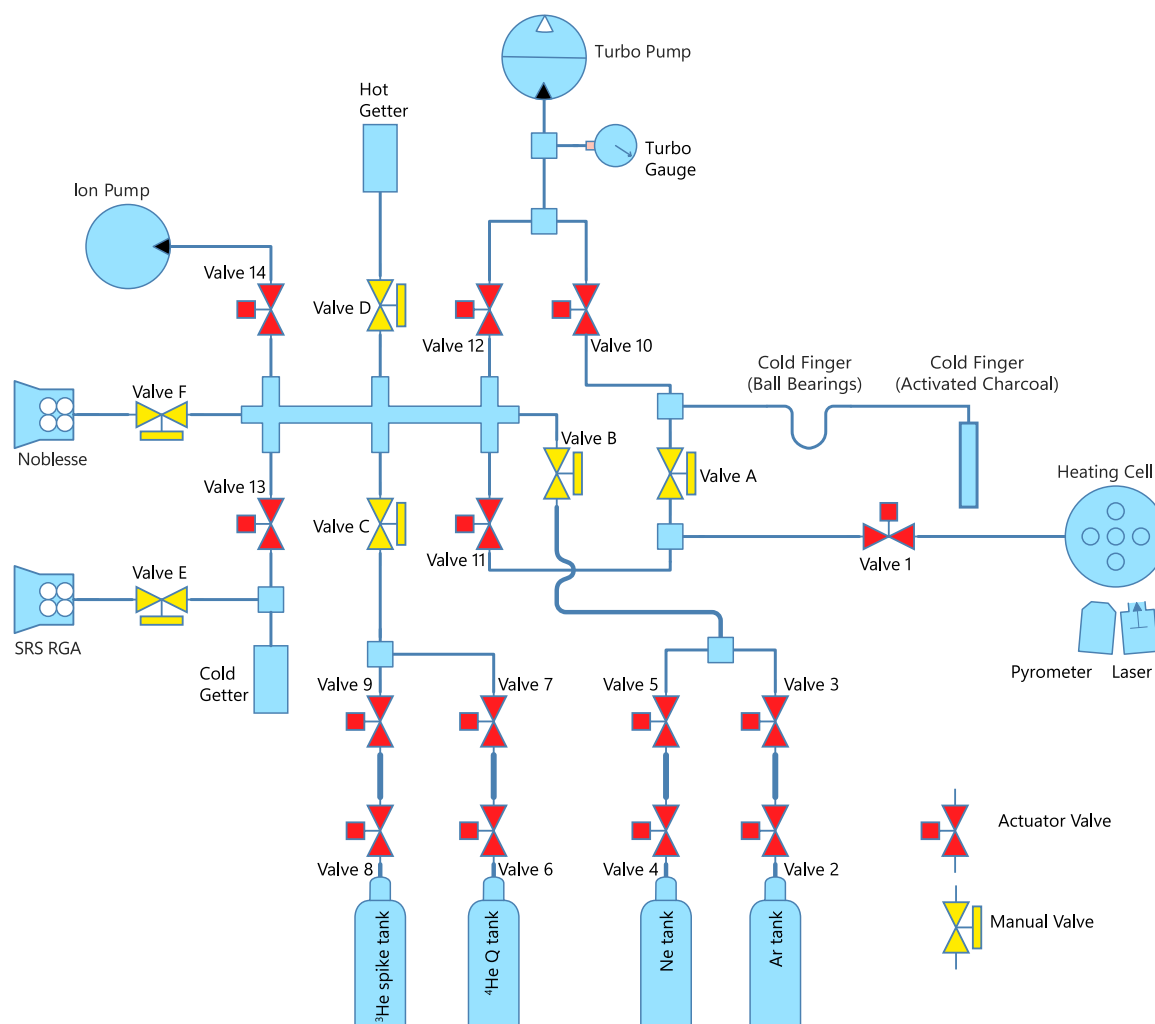


Figure 1: The UCL Oxide Nobel Gass Extraction Line

It consists of a laser/pyrometer system use to heat samples and a mass spectrometer used to measure the resultant gasses. The extraction line includes a ^3He spike tank, a ^4He standard tank as well as an Ar and Ne tank. There is a cryotrap which consists of two cold fingers that can be added into the line as well.

Valve	Description
Valve A	Cryotrap Isolator
Valve B	Ar / Ne Tank Isolator
Valve C	Helium Tanks Isolator
Valve D	Hot Getter Isolator
Valve E	SRS Mass Spectrometer Isolator
Valve F	Noblesse Mass Spectrometer Isolator

Table 1: Manual Valves

Valve	Description
Valve 1	heating cell
Valve 2	Ar tank pipette input
Valve 3	Ar tank pipette output
Valve 4	Ne tank pipette input
Valve 5	Ne tank pipette output
Valve 6	4He Q tank pipette input
Valve 7	4He Q tank pipette output
Valve 8	3He spike tank pipette input
Valve 9	3He spike tank pipette output
Valve 10	turbo to cryotrap
Valve 11	input to manifold
Valve 12	turbo to manifold
Valve 13	SRS RGA
Valve 14	Ion Pump
Valve 15	spare

Table 2: Automated valves

Valve Controller Description

The valve and laser control software is a Python (Python Software Foundation 2023) application that manages up to 15 solenoid valves. The Raspberry Pi 4B (Raspberry PI Foundation 2020) controls these valves using a 40-way connector for the necessary GPIO pins and driver boards. Each valve operates via a dedicated GPIO line. Since the valves require 24V DC and the Pi uses 3.3V signals, a driver circuit steps up the voltage and current.

The system is managed through a RESTful API on port 80, accepting valid JSON messages. A 128 character API-Key is generated when the application is first run to secure the API.

Safety logic prevents simultaneous opening of pipette input and output valves (valve pairs 2-3, 4-5, 6-7). A single command can close all valves in case of an issue.

Additionally, the controller monitors two vacuum gauges, a Pfeiffer MPT200 gauge which is connected to a RS485 USB to serial adapter and a Gamma Vacuum SPc ion pump using a CH340 USB to RS232 controller. Both gauges are polled every 4 seconds to read the gas pressure.

Input Messages

Json messages in the following formats are accepted:

Open a single valve (nn):

```
{  
  "item": "valvenn",  
  "command": "open"  
}
```

Close a single valve (nn):

```
{  
  "item": "valvenn",  
  "command": "close"  
}
```

Close all valves:

```
{  
  "item": "closeallvalves"  
  "command": ""}
```

Get Valve Status: (for all valves)

```
{  
  "item": "valvestatus",  
  "command": ""  
}
```

Read Pressures:

```
{  
  "item": "getpressures",  
  "command": "read"  
}
```


Output Messages

Following a “valvestatus” command any valid command, the following data is returned: (ss = open or closed, xx = on or off)

```
{
  "status": "closed",
  "valve": 1
},
{
  "status": "closed",
  "valve": 2
},
{
  "status": "closed",
  "valve": 3
},
{
  "status": "closed",
  "valve": 4
},
{
  "status": "closed",
  "valve": 5
},
{
  "status": "open",
  "valve": 6
},
{
  "status": "open",
  "valve": 7
},
{
  "status": "open",
  "valve": 8
},
{
  "status": "open",
  "valve": 9
},
{
  "status": "open",
  "valve": 10
},
{
  "status": "closed",
  "valve": 11
},
{
  "status": "closed",
  "valve": 12
},
{
  "status": "closed",
  "valve": 13
},
{
  "status": "closed",
  "valve": 14
},
{
  "status": "closed",
  "valve": 15
}
```

Following a “getpressures” command, the following data is returned:

```
{  
  "pressure": 4.17e-08,  
  "pump": "turbo"  
},  
{  
  "pressure": 1.4e-09,  
  "pump": "ion"  
}
```

Pressures will be in floating point notation when connected to the gauge. When disconnected they will return zeros

Valve to GPIO assignment

Valve	Designation	Connector		Designation	Valve
	3v3	1	2	5v	
	GPIO 02	3	4	5v	
	GPIO 03	5	6	GND	
	GPIO 04	7	8	GPIO 14	
	GND	9	10	GPIO 15	
Valve 5	GPIO 17	11	12	GPIO 18	Valve 4
Valve 3	GPIO 27	13	14	GND	
Valve 2	GPIO 22	15	16	GPIO 23	Valve 1
	3v3	17	18	GPIO 24	Valve 10
	GPIO 10	19	20	GND	
Valve 9	GPIO 09	21	22	GPIO 25	
Valve 8	GPIO 11	23	24	GPIO 08	
	GND	25	26	GPIO 07	
	GPIO 00	27	28	GPIO 01	
	GPIO 05	29	30	GND	
	GPIO 06	31	32	GPIO 12	Valve 7
Valve 6	GPIO 13	33	34	GND	
Valve 15	GPIO 19	35	36	GPIO 16	Valve 14
Valve 13	GPIO 26	37	38	GPIO 20	Valve 12
	GND	39	40	GPIO 21	Valve 11

Table 3: GPIO to valve assignments

Table 2 shows GPIO channels on a Raspberry Pi 4B, channels in green are available and remain at 0v during the Pi boot up sequence, it is important that no lines are used that could cause the laser to turn on or a valve to open before the software is ready.

Driver circuits and power supply

The Raspberry Pi 4 operates with a 5V power supply delivered through a USB-C connector. A power supply rated at 3A is necessary for optimal performance. The internal operating voltage of the Raspberry Pi and the output voltage of the GPIO connectors are both 3.3V.

The valves are activated by supplying 24V DC to the solenoid. Upon removal of the power, the valves will close due to a spring mechanism. It is assumed that up to 10 valves may be actuated simultaneously during a pump-down cycle; therefore, the power supply must be capable of delivering at least 1.0A continuously. To ensure reliable provision of this current, a 24V power supply with a continuous rating of 3.0A is required.

To safeguard the Raspberry Pi from potential voltage spikes generated by solenoids during valve closure or due to a malfunctioning component, the use of an opto-isolator is recommended. This will ensure complete electrical isolation between the Raspberry Pi and the driver circuit.

The solenoids on the valves will be operated using Metal Oxide Field Effect Transistors (MOSFETs). MOSFETs exhibit a very low forward voltage once activated and possess a high power capacity (Inchange Semiconductor 2016).

Since the valves are operated by inductive solenoids, removing power from the solenoids can potentially induce voltage spikes into the circuit. Therefore, a Schottky diode is used to protect the MOSFET and opto-isolator from any potential damage.

Each valve driver board has 5 channels, so the valve controller requires 3 driver boards to control the 15 valves (Figure 2).

The driver circuit must be enclosed in a housing that includes external connectors for the mains supply, Ethernet cable, and 12 valve cables. Ventilation holes are necessary to prevent the Raspberry Pi, power supplies, or driver MOS-FET transistors from overheating. Additionally, indicator lights should display the status of the 5V PSU for the Raspberry Pi and the 24V valve PSU.

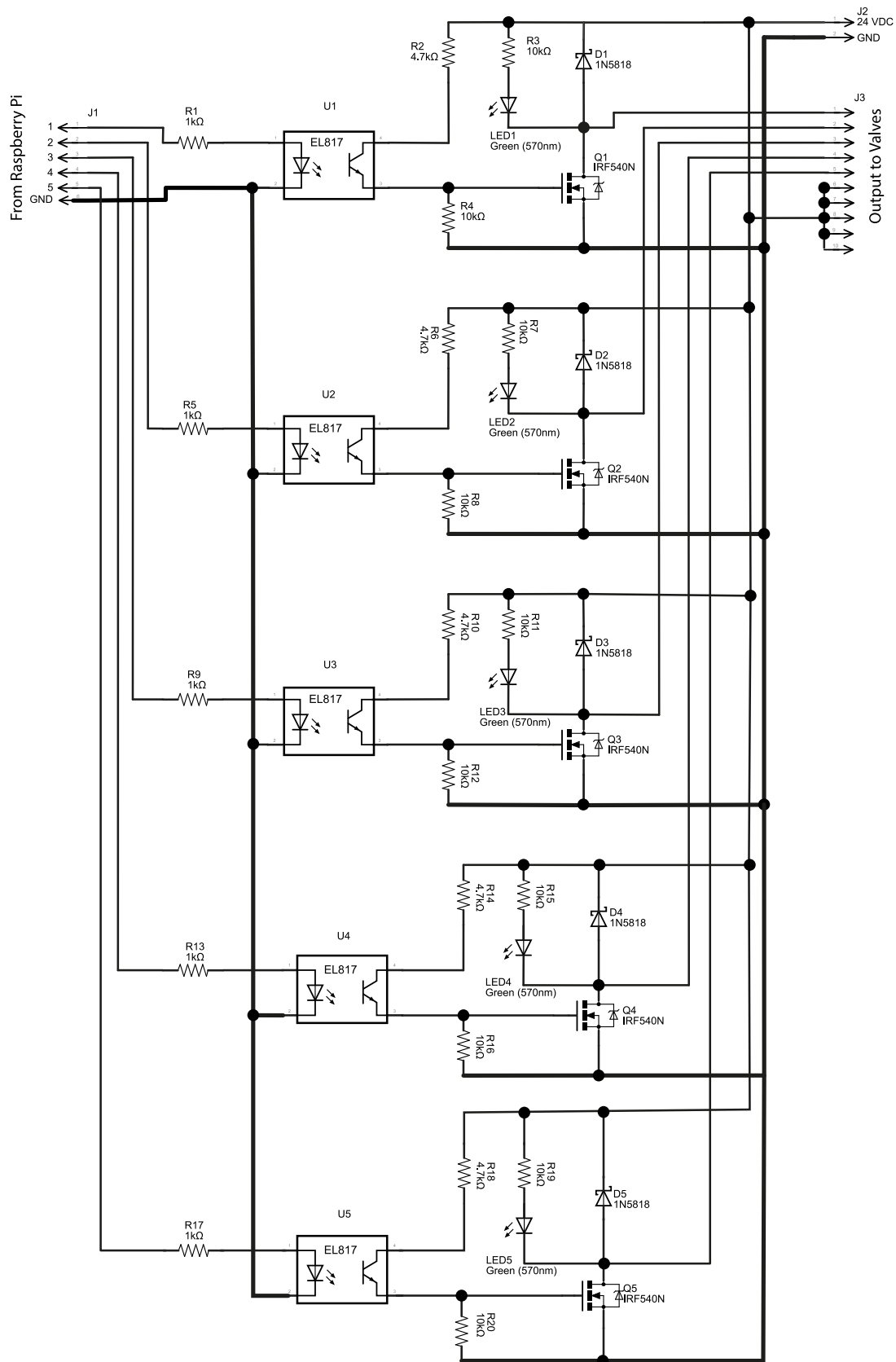


Figure 2: Schematic for a valve driver board

Web interface

In addition to accessing the status of the valves via the RESTful API, a read-only web interface will be available. This interface can be accessed directly from a browser to view the valve status, as well as application and web server logs.

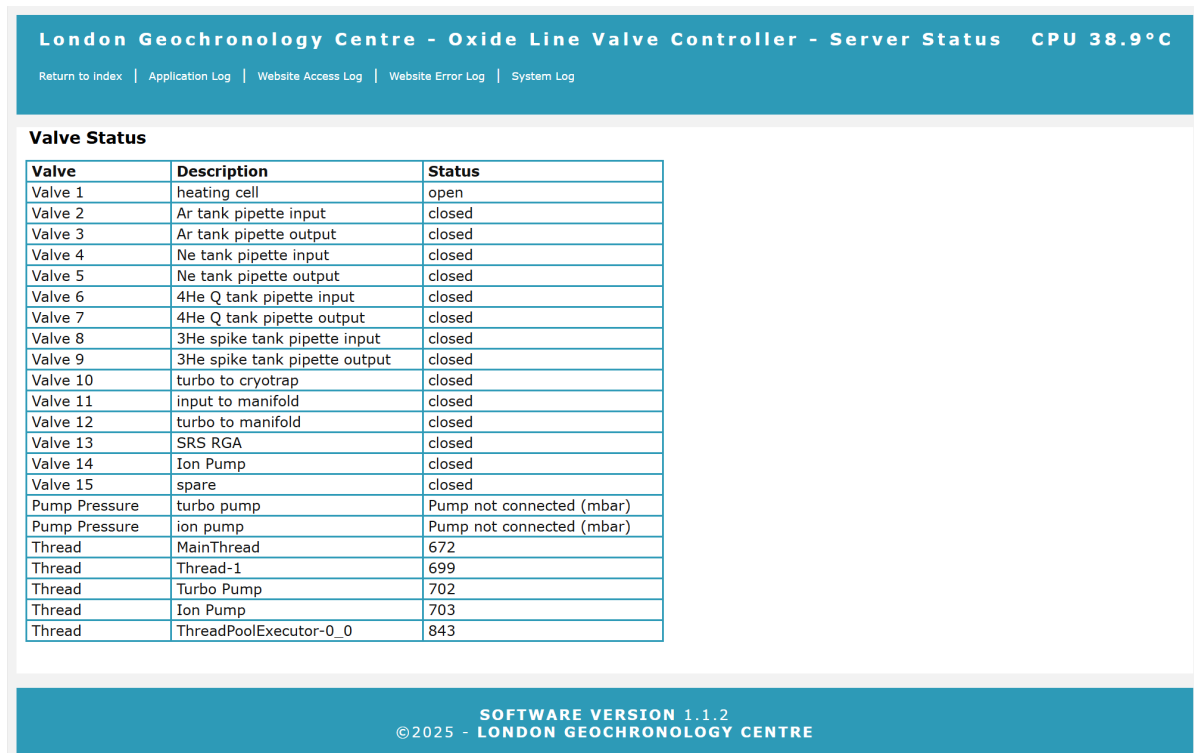


Figure 3: Web status page

Interfaces

9 Pin Valve Cable Specification

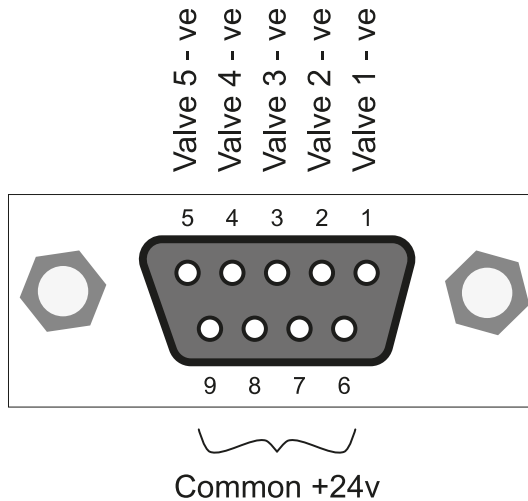


Figure 4: Valve controller 9 pin D plug assignments

Three 9 pin D plugs connected to 5 individual cables each, terminated with an SMC 24v latching power plug. 3 plugs will be required for the 15 valves.

9 Pin RS232 Cable specification

The connection from RS232 port a standard 9 pin D Socket, a null modem switching pins 2 and 3 (receive and transmit) is used to connect to the Gamma Vacuum SPc gauge..

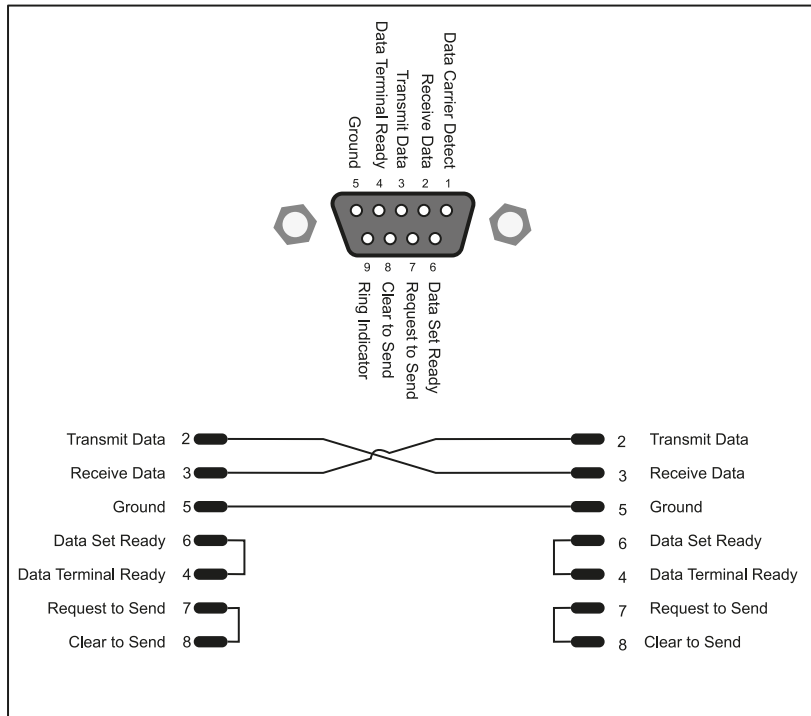


Figure 5: RS232 D Plug Assignment and null modem cable

4 Pin RS485 Cable Specification

The cable intercepts the RS485 signal sent from the Pfeiffer MT200 gauge to the Pfeiffer controller. The Pfeiffer equipment uses an M12 4 pin instrumentation connector.

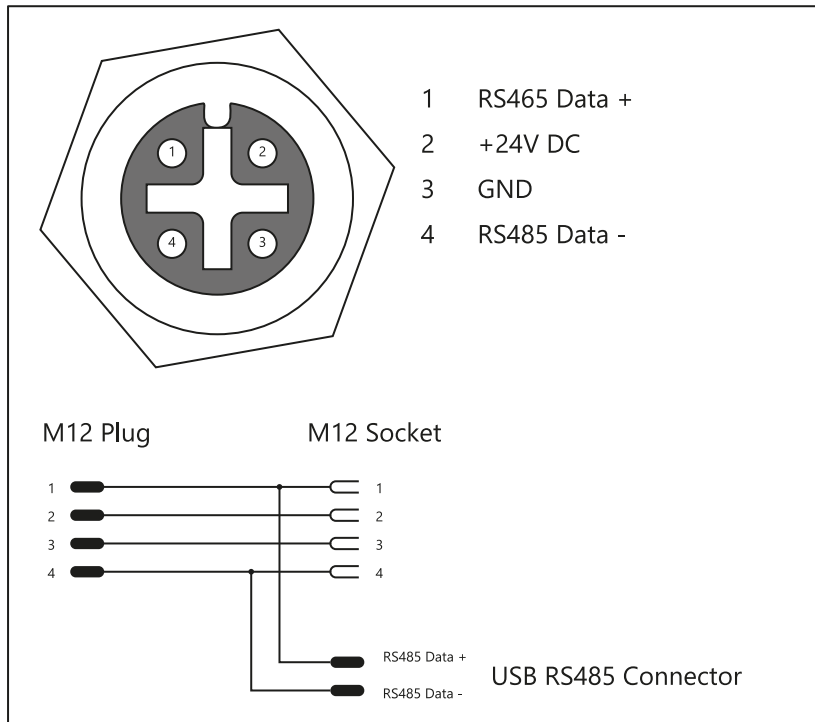


Figure 6: RS485 M12 connector layout and intercept cable

USB Connections

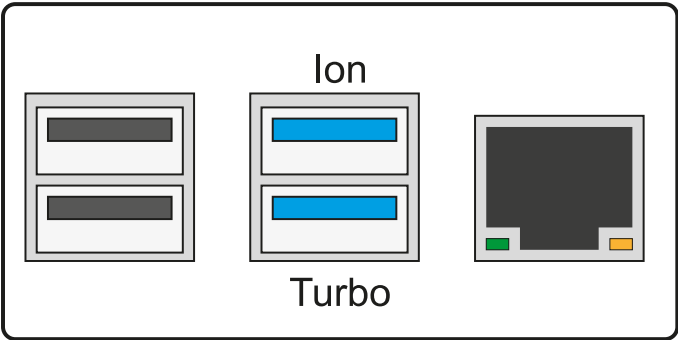


Figure 7: USB connections to the Valve Controller Raspberry Pi

Connection	USB	Serial Device	Equipment Connector
Ion Pump	USB3 port 0	/dev/ttyUSB0	RS232 Socket (female)
Turbo Pump	USB3 port 1	/dev/ttyUSB1	RS232 Plug (Male)

Table 4: External Equipment Connections

Gauge Data Request Commands

Equipment	First Command (Hexadecimal)	Second Command (Hexadecimal)
Pfeiffer Vacuum TPG 1 (Tank and Turbo)	50 52 31 0D (Pressure Gauge 1)	05 (Enquiry)
Gamma Vacuum SPc (Ion Pump)	7E 20 32 35 20 30 42 20 30 30 0D (read pressure from ID 05)	

Table 5: Commands to Vacuum Gauges

Raspberry Pi Operating System Installation

Operating system

Use BalenaEtcher to install the latest version of the bookworm-lite 64bit operating system onto a 32Gb MicroMMC card.

Run the `sudo raspi-config` command to:

- ♦ enable ssh
- ♦ disable Serial
- ♦ disable 1²C bus
- ♦ set the password on the pi user account to something other than “raspberry”
- ♦ change the GPU memory to 16Gb

Raspberry Pi Software installation

Run Operating System Updates

Run `sudo apt update`

Run `sudo apt upgrade`

Install git and gh

Run `sudo apt install git`

Run `sudo apt install gh`

Run `gh auth login`

- What account do you want to log into? GitHub.com
- What is your preferred protocol for Git operations? HTTPS
- Authenticate Git with your GitHub credentials? (Y/n) y
- How would you like to authenticate GitHub CLI? Paste an authentication token

Configure git and download the valve controller code

Create a GitHub folder

Run `mkdir github`

Run `cd github`

Clone the repo

Run `git clone https://github.com/westerlymerlin/UCL-Oxide-Valve-Controller.git`

Copy the files to the home folder

Run `cp -r ~/github/UCL-Oxide-Valve-Controller/* ~/`

Set the execute flag on the scripts

Run `chmod 755 ~/bin/*`

Copy the Raspberry pi config files to the etc folder

Run `sudo cp -r ~/raspberry-pi/etc/* /etc`

Reboot the Raspberry pi

Run `sudo reboot`

After the reboot a warning banner will appear at ssh logon.

Create a python venv

Run `python3 -m venv .venv`

Change to the venv

Run `source activate`

Install the Required Python Libraries

Run `pip install -r requirements.txt`

Gunicorn for Python service

Run `sudo systemctl enable gunicorn`

Run `sudo systemctl start gunicorn`

Nginx installation

Run `sudo apt install nginx`

Change directory to `/etc/nginx/sites-enabled/`

Run `sudo rm default`

Run `sudo ln -s /etc/nginx/sites-available/valvedata`

Reboot the Raspberry pi

Run `sudo reboot`

If flask is installed, the python files are in the `/home/pi` directory, gunicorn3 is installed and configured and nginx is installed and configured the web service should be running, and the site will be accessible on `http://ip address of the server`

References

- Inchange Semiconductor (2016) IRF540N. Available online: <http://www.iscsemi.com> [Accessed].
- Python Software Foundation (2023) *Python 3 Programming Language*. Online. Available online: <https://www.python.org> [Accessed March 2023].
- Raspberry PI Foundation (2020) *Raspberry PI Model 4B Reference*. Available online: <https://www.raspberrypi.org/products/raspberry-pi-4-model-b/> [Accessed July 2020].